

Full paper / Mémoire

Spatially resolved operando infrared analysis of a microstructured catalytic surface for CO oxidation over Pt based catalysts

Cécile Daniel ^a, Marie-Olive Clarté ^a, Hélène Provendier ^a, André C. Van Veen ^a,
Yves Schuurman ^a, Bruno J. Beccard ^b, Claude Mirodatos ^{a,*}

^a IRCELYON, Institut de recherches sur la catalyse et l'environnement de Lyon, Université Lyon 1, CNRS, UMR 5256,
2, avenue Albert-Einstein, 69626 Villeurbanne, France

^b Thermo Fisher Scientific, 16, avenue du Québec, SILIC 765, 91963 Courtaboeuf cedex, France

Received 15 June 2008; accepted after revision 4 September 2008

Available online 28 February 2009

Abstract

A new experimental reactor concept was developed for operando and spatially resolved DRIFTS surface analysis over a microstructured catalytic reactor. Both gas and surface adspecies concentration profiles can be established along the channels coated with catalyst.

The oxidation of carbon monoxide was studied from 25 to 300 °C over a Pt/CeO₂–Al₂O₃ catalyst. For these operating conditions nearly isothermal conditions were observed for the microstructured reactor, with a continuous increase in conversion along the microchannels. In contrast, an important thermal effect with characteristic light-off was obtained in the fixed bed reactor, underlining the unique behaviour of the microstructured system.

Modelling of the CO conversion and surface adspecies concentration as a function of temperature along the channels was performed on the basis of a simple reaction mechanism. Kinetic parameters depending on the catalyst composition were obtained, in good agreement with literature data. **To cite this article:** C. Daniel et al., *C. R. Chimie* 12 (2009).

© 2009 Académie des sciences. Published by Elsevier Masson SAS. All rights reserved.

Résumé

Un nouveau type de réacteur a été conçu pour une analyse DRIFTS operando et résolue spatialement d'une surface catalytique structurée. Les concentrations en gaz et en espèces adsorbées peuvent être ainsi mesurées le long de microcanaux revêtus d'une couche catalytique. Un réacteur conventionnel DRIFTS en lit fixe a été testé dans les mêmes conditions pour comparaison.

La réaction d'oxydation du monoxyde de carbone a été étudiée de 25 à 300 °C sur un catalyseur Pt/CeO₂–Al₂O₃. Un comportement quasi isotherme a été observé pour ce réacteur microstructuré, avec une augmentation régulière de la conversion le long des canaux. A l'inverse un effet thermique marqué avec un allumage rapide de la réaction a été observé dans le cas du réacteur à lit fixe, confirmant les excellentes propriétés de transfert thermique du réacteur microstructuré.

* Correspondent Author.

E-mail address : claude.mirodatos@ircelyon.univ-lyon1.fr (C. Mirodatos).

Une modélisation de la conversion du CO et des concentrations en espèces adsorbées le long des micro-canaux a été réalisée sur la base d'un modèle cinétique simplifié. Les paramètres cinétiques qui dépendent de la composition du catalyseur se sont avérés en bon accord avec les données de la littérature. **Pour citer cet article : C. Daniel et al., C. R. Chimie 12 (2009).**

© 2009 Académie des sciences. Published by Elsevier Masson SAS. All rights reserved.

Keywords: Operando DRIFTS spectroscopy; CO oxidation; DRIFTS transparent microreactor; Pt/CeO₂-Al₂O₃

Mots-clés : Spectroscopie infrarouge à réflexion diffuse en mode operando ; oxydation du CO ; réacteur microstructuré transparent à l'infrarouge ; Pt/CeO₂-Al₂O₃

1. Introduction

Within the key perspective of process intensification, microstructured reactors appear as a highly promising technology [1] and perspectives of mass production are announced by technology providers [2,3]. These reactors are generally based on the assembly/stacking of microstructured metallic platelets coated with suitable catalysts. Due to their high heat and/or mass transfer, low pressure drop and good phase contacting, they are particularly adapted to the large domain of catalytic oxidation reactions. Specificities for catalyst wash coating, platelet assembly and operating conditions have been carefully considered for implementing correctly these new structured reactors [4]. However, a major challenge remains to investigate the catalytic process along the microchannels in view of further modelling and better control of the operating conditions. For that purpose, an access to effective gas and surface adspecies concentration profiles is required, rather than establishing it only by calculation. Until now, only some attempts have been reported to establish a gas concentration profile by sampling the effluents with a micro-capillary along structured reactors (in general monolith based), completed with a thermal mapping for detecting hot spots along the catalytic bed [5].

In this work, a new experimental reactor concept is described, which allows true surface analysis by Diffuse Reflectance Infrared Fourier Transform Spectroscopy (DRIFTS) measurements over a microstructured reactor, under reaction conditions through an adapted window without affecting the flow along the micro-channels. The reaction chosen to validate this concept is the oxidation of carbon monoxide, which represents a key step in the process of hydrogen purification before feeding PEM fuel. A Pt/CeO₂-Al₂O₃ catalyst, known to be active for that reaction, was selected for this preliminary operando DRIFTS analysis.

In order to evidence specificities of the micro-structured system, the related results were compared to those obtained in a conventional fixed bed DRIFTS reactor with a similar catalyst in powder form and under closed operating conditions.

2. Experimental

2.1. DRIFTS reactors

2.1.1. Optically accessible DRIFTS microstructured reactor

Fig. 1A shows a schematic design of the optically accessible microstructured reactor. The catalyst is

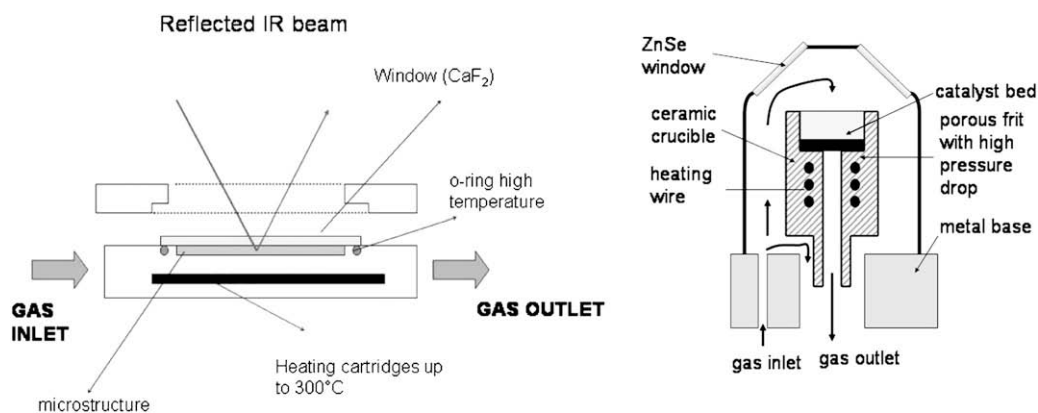


Fig. 1. Design of the optically accessible DRIFTS reactors. (A): microstructured (in house), (B): fixed bed (from Spectratech[®]).

coated as a thin layer onto a microstructured metal platelet. This platelet is placed into a reactor housing designed to be adapted to a xy table of an Infrared optical microscope “Nicolet Continuum” linked to an infrared spectrometer “Nicolet Nexus” equipped with an MCT detector [6]. This set-up is connected to a conventional gas supply ensemble and an analysis set-up (micro-gas chromatograph Agilent 3000 instrument) equipped with a channel module used for H₂, O₂, N₂, CH₄ and CO measurements (MolSieve 5A, BF PPU Plot, with Ar as carrier gas) and a channel for CO₂ and H₂O (PPU, BF PPQ, with He as carrier gas). The housing can be heated up to 300 °C by electrical cartridges. Its design allows a laminar flow well distributed among the parallel channels, as demonstrated by DFT calculations (not reported here). The CaF₂ window is directly in contact with the platelet, avoiding significant by-pass between the channels. By moving automatically the reactor under the infrared beam, the spatially resolved DRIFTS analysis is carried out along the reactor channels under reaction conditions. For this reactor, the optical length of the IR beam within the reacting zone, i.e. the geometrical distance between the reactor window and the catalyst adlayer, is evaluated to be about 0.5 mm (see Fig. 1). The IR signal was collected in the form of absorbance ($\log(1/R)$ transform) which is the most adapted to the present case of poorly absorbing adsorbates to obtain DRIFTS data proportional to the adsorbate surface concentrations [7]. Background single beam spectra were acquired by focussing on a catalyst free zone of the microreactor to account for all internal and external non-catalytic contributions. Though a minor specular reflection contribution cannot be discarded, it is assumed that the signal comes essentially from diffuse reflectance.

2.1.2. Fixed bed DRIFTS reactor

The DRIFTS analysis carried out at the surface of a fixed bed reactor loaded with the same catalyst but packed under powder form was performed by using an *in situ* DRIFTS cell from Spectratech [8] (Fig. 1B).

About 20 mg of powder catalyst sieved at 200–400 μm was loaded in the cell on a layer of silicon carbide to improve the thermal exchange within the fixed bed. A CO/O₂/Ar mixture was passed over the catalyst at different temperatures.

For this reactor, the optical length of the IR beam within the reacting zone is evaluated to be about 10 mm. A partial by-pass of the fixed bed as indicated in Fig. 1 prevented full conversion to occur, as already

stated by Meunier et al. [9]. However, the effective conversion levels were measured in a conventional tubular fixed bed reactor under the same operating conditions. The same mode of spectra acquisition was used as for the microstructured reactor.

2.2. Catalysts preparation and platelets coating for the microstructured DRIFTS reactor

The structured platelet used in this study was made of Aluchrom[®] (Iron–chromium–aluminium resistance alloy), cleaned and annealed at 1000 °C before coating. An alumina sublayer was deposited first by an optimized wash coating recipe (using a solution of acrylic acid mixed with a suspension of alumina, acetic acid and boehmite), followed by drying at room temperature and slow calcination (1 K/min) up to 1000 °C. Ceria was then impregnated from a solution of Ce(NO₃)₃ and the mixed oxide layer slowly calcined up to 550 °C. About 2 wt.% of platinum was then added by impregnation from a solution of Pt(NH₃)₄(OH)₂ followed by a slow calcination (1 K/min) up to 550 °C [10].

2.3. Catalysts preparation for the fixed bed DRIFTS reactor

The powder catalyst used in the fixed bed DRIFTS reactor was prepared by a sol–gel method with a formulation close to the one used for the structured microreactor [11,14].

2.4. DRIFTS kinetic measurements and gas analysis of effluents

For both DRIFTS reactors (fixed bed or microstructured), the catalysts were reduced under a flow of 50% H₂ in He from room temperature to 270 °C at 10 °C min⁻¹ and kept 2 h at 270 °C before cooling under He.

In a typical catalytic test for investigating temperature effect, the reactors were heated from 25 °C to 300 °C, and eventually ramped down to room temperature. The DRIFTS spectra were recorded continuously and the bands' intensity calculated after proper background treatments.

2.4.1. Microstructured DRIFTS reactor

During catalytic testing for mapping gas and adspecies concentrations along the microstructured channel, the operating conditions (flow rates and temperature) were maintained constant and the microstructured reactor displaced regularly under the IR beam until the whole channel was investigated

(6 analyses per channel). The gas inlet composition was fixed at $\text{CO}:\text{O}_2:\text{He} = 50\%:25\%:25\%$. This high concentration of reacting species was chosen to maximize the DRIFTS signals due to the short optical length (around 0.5 mm).

2.4.2. Fixed bed DRIFTS reactor

This reactor offering an optical length about 20 times larger than for the microstructured one, the gas inlet composition was fixed as either $\text{CO}:\text{O}_2:\text{Ar} = 0.5\%:0.6\%:98.9\%$, to maximize the absorbance of adsorbed species and minimize the one of gas phase reactants and products, or $\text{CO}:\text{O}_2:\text{He} = 50\%:25\%:25\%$ to compare strictly with the conditions selected for the microstructured system.

During all these operations and for both DRIFTS reactors, the outlet gas concentration was analysed by the on-line micro-gaseous chromatograph (similar to the microstructured system).

3. Results

3.1. CO oxidation in the fixed bed DRIFTS reactor

The spectra obtained during the CO oxidation reaction in the fixed bed DRIFTS reactor on the $\text{Pt}/\text{CeO}_2\text{--Al}_2\text{O}_3$ catalyst were acquired in the temperature range 20–300 °C by ramping up the system at 5 °C/min, with 8 min at each constant temperature interval (20, 50, 100, 150, 200, 250, 300 °C). They revealed two domains of temperature corresponding to two distinct states of the catalytic surface, as shown in Fig. 2.

- (i) from room temperature to about 180 °C, the band of CO linearly adsorbed on Pt phase at 2069 cm^{-1} (named L-CO in Fig. 2) remains rather narrow and stable, as well as the doublet of gaseous CO at

$2100\text{--}2175\text{ cm}^{-1}$ (not reported in Fig. 2, due to its small intensity). In the mean time the doublet at $2360\text{--}2330\text{ cm}^{-1}$, characteristic of gaseous CO_2 (named CO_2 gas in Fig. 2), develops very slowly from 100 °C to above. Other adspecies characteristic of the basic ceria–alumina support (formates and carbonates adspecies) were also observed in this temperature range within the range $1622\text{--}1554\text{ cm}^{-1}$.

- (ii) from 180 to 300 °C, the doublet of gaseous CO and the linearly adsorbed CO band at 2065 cm^{-1} are no longer detected, while the gaseous CO_2 doublet is now large and stable with temperature. This abrupt change in surface and gas concentration obviously reflects the light-off and run away of the CO oxidation, leading to the full conversion of CO into CO_2 (at least for the flow crossing the reactor), as checked in a conventional fixed bed reactor without by-passing effect. No major changes were observed for the support adspecies.

This temperature light-off (leading to a higher temperature of the catalyst bed as compared to the reactor temperature) is confirmed by observing a marked hysteresis between ramping up and ramping down of temperature (A). This latter effect is the characteristic of a quasi autothermal system where the heat produced by the reaction overcompensates the heat losses by convection and/or conduction.

3.2. CO oxidation in the microstructured DRIFTS reactor

Fig. 3 illustrates the DRIFTS mapping along a central microchannel (every 6 mm) under reacting conditions (here at 300 °C). Due to the small optical length (0.5 versus 10 mm) and the small amount of

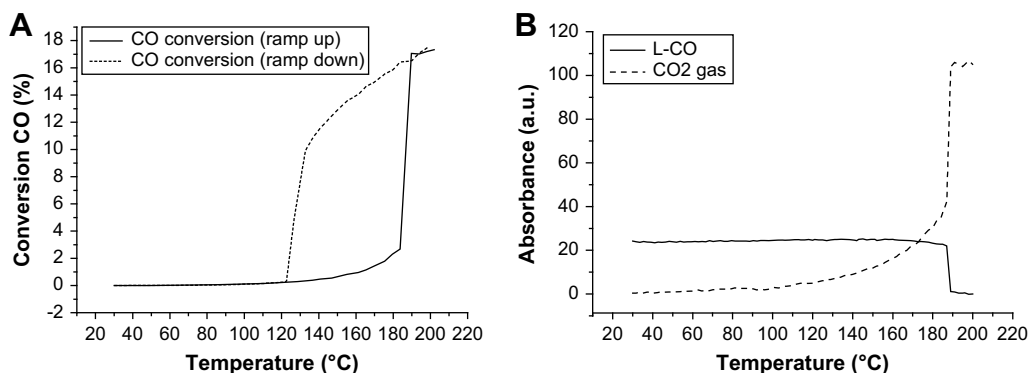


Fig. 2. Changes in CO conversion at reactor outlet (A) and in CO_2 gas phase and linear carbonyls concentration (B) versus temperature for the $\text{Pt}/\text{CeO}_2\text{--Al}_2\text{O}_3$ in the fixed bed DRIFTS reactor. $\text{CO}:\text{O}_2:\text{He} = 50:25:25$.

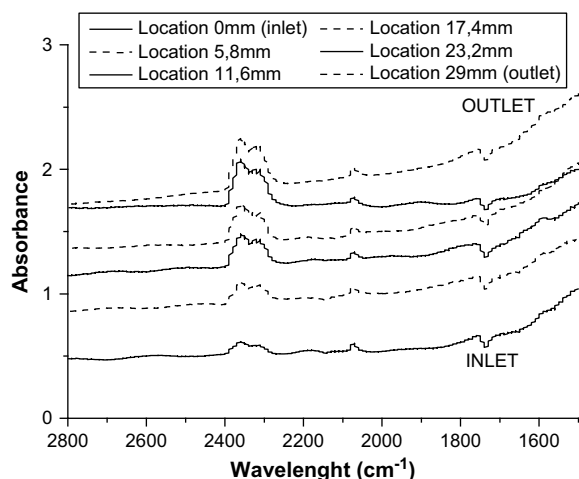


Fig. 3. Changes in DRIFTS spectra from inlet to outlet in the microstructured reactor during the oxidation of carbon monoxide at 300 °C over Pt/CeO₂–Al₂O₃.

reflecting material in that device, only the most pronounced bands (L-CO, gaseous CO and CO₂ and adsorbed formates/carbonates lumped together) are quantified.

Fig. 4 reports the changes in gaseous CO₂ band intensity along the reactor axis for the Pt/CeO₂–Al₂O₃ catalyst. This concentration is expressed as a local CO₂ flow rate after calibration done at the reactor outlet where the effective CO₂ flow rate (obtained by mass chromatography) is compared to the DRIFT response of the CO₂ gas IR band.

Quasi linear increases of CO₂ concentration with temperature are observed, whatever the temperature.

As seen in Fig. 5, in contrast with the previous linear increases in CO₂ gas concentration along the microchannels, no significant change in surface concentration was observed for the main detected

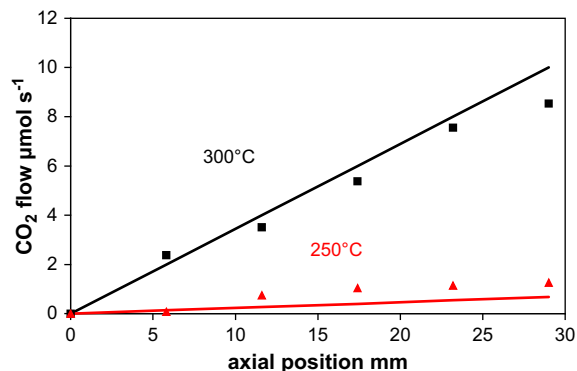


Fig. 4. Experimental (symbols) and calculated (lines) CO₂ concentration (expressed as local flow rate) profiles along a microchannel at 250 and 300 °C over the Pt/CeO₂–Al₂O₃ catalyst.

adspecies, within the uncertainty of spectrum acquisition. Note here that the latter was mainly due to: i) the weakness of the IR signal coming from the limited amount of reflecting material (see Fig. 3); and ii) the lack of homogeneity of the coated adlayer (thickness, conformity, microcracks, surface structural defects, etc), revealed by the optical microscope and quite in line with the spot resolution (50 × 50 μm). Thus, formates and carbonates adspecies were lumped together within the range 1622–1554 cm⁻¹ and linear carbonyl adsorbed on platinum particles were detected within the range 2084–2053 cm⁻¹. The above statement indicated a constant coverage of the catalyst all along the channels, within the explored temperature range (225–300 °C).

3.3. Changes in conversion at reactor outlet

Fig. 6 reports the changes in CO conversion with temperature measured at the outlet of the DRIFTS microstructured reactor. At variance with the CO conversion curves obtained with the fixed bed DRIFTS reactor, no light-off effects (abrupt increase of conversion around 180 °C and hysteresis loop between ramp up and down temperature changes) were observed with the microstructured reactor: the conversion increased in a monotonic way, with the same values when ramping the temperature up and down, as expected from an isothermal and kinetically controlled process.

4. Discussion

4.1. Comparison between the microstructured and fixed bed DRIFTS reactors

As depicted in Figs. 2 and 6, thermal effects were observed only in the fixed bed DRIFTS reactor. This is a characteristic of a reaction light-off and run away (generally ascribed to heat transfer limitation occurring in non-isothermal systems), leading to full conversion of linear carbonyl adspecies on Pt particles, with the simultaneous disappearance of the gas phase CO and the abrupt increase in gaseous CO₂ (Fig. 2). The two temperature domains observed experimentally correspond respectively to two distinct states of the active surface: i) below 180 °C, full coverage of Pt particles by linear carbonyls, strongly inhibiting the oxygen activation; and ii) above 180 °C, the Pt particles are mainly covered with activated oxygen, immediately reacting with CO as soon as colliding with the Pt surface.

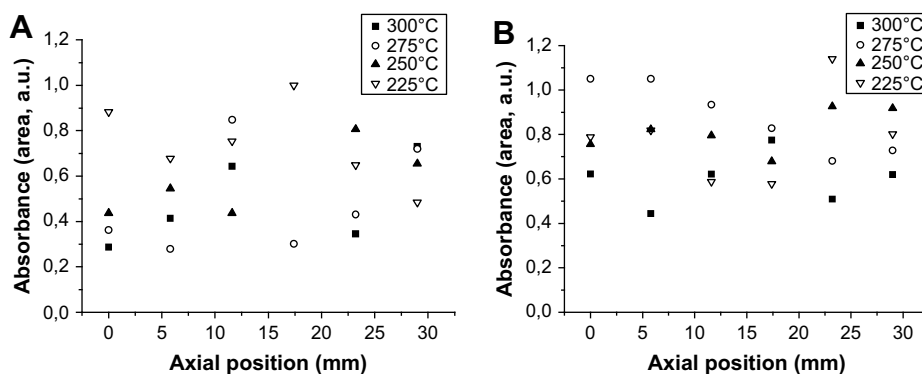


Fig. 5. Changes in (A) formates and carbonates adspecies (lumped within $1622\text{--}1554\text{ cm}^{-1}$) and (B) linear carbonyl on platinum ($2084\text{--}2053\text{ cm}^{-1}$) along the microstructured channels between 225 and 300 °C for the Pt/CeO₂–Al₂O₃ catalyst.

In contrast, over the whole 25–300 °C temperature range, the microstructured reactor leads to a monotonous change in CO conversion with temperature (Fig. 6). This demonstrates the excellent heat transfer in this microstructured system which prevents any run away of the CO oxidation, confirming this unique property of microstructured reactors [1,4].

4.2. Specificity of the Pt/CeO₂–Al₂O₃ catalyst in the microstructured reactor

As demonstrated by the DRIFTS analysis carried out in the microstructured reactor (illustrated in Fig. 5), the coverage of platinum particles by linear carbonyls remains high all over the explored reacting conditions (little changes between room temperature and 300 °C), due to the efficient heat transfer with the metallic substrate, evacuating continuously the heat produced by the CO oxidation. Therefore, it can be inferred that under these conditions, no Pt sites are available for oxygen activation, which means, as generally agreed in the literature [13], that the main source of active oxygen is ensured by ceria, able to activate gaseous oxygen for feeding continuously the Pt/ceria interface where the CO oxidation can proceed. Note that the proper intrinsic activity of the ceria/alumina support has been found negligible from specific experiments not reported here.

4.3. Modelling of CO oxidation over Pt/CeO₂/Al₂O₃ in the microstructured reactor

Campman [12] and Nibbelke et al. [13] studied extensively CO oxidation over a three-way Pt/CeO₂/Al₂O₃ catalyst. They proposed first the following mechanism on the Pt phase:



Assuming that the O₂ adsorption is rate-limiting and that the CO adsorption is equilibrated give the following rate equation:

$$r_{\text{CO}} = \frac{2L_t k_2^f P_{\text{O}_2}}{(1 + K_1 P_{\text{CO}})} \text{ mol/s/kg}_{\text{cat}} \quad (5)$$

The following a priori values of the kinetic parameters for Pt were selected from literature.

The priori values of the parameters for Pt, mainly from surface science studies (Table 1), were substituted in rate Equation (5) and incorporated into a one-dimensional plug flow reactor model. The corresponding differential equation was integrated

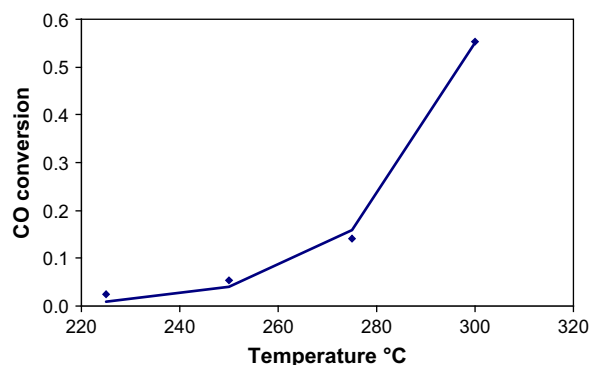



Fig. 6. Changes in CO conversion with temperature at the outlet of the DRIFTS microstructured reactor for Pt/CeO₂–Al₂O₃ catalyst. Experimental (symbols) and modelled (continuous line) data.

Table 1
A priori values of the kinetic parameters for Pt.

Parameter	Value	References
L_t	0.04 mol/kg _{cat} (~ 60% dispersion)	Germani et al. [14]
k_2^f	$0.2\text{--}2 \times 10^7/\text{bar s}$ ($s^0 \approx 0.01\text{--}0.1$)	Luntz [15]
E_2^f	0 kJ/mol	Herz and Marin [16]
K_1	$2 \times 10^{-6}/\text{bar s}$	Zhdanov et al. [17]
E_1	-130 kJ/mol	Herz and Marin [16]

numerically over the reactor length and the value of k_2^f was adjusted within the range of values given in Table 1 to give an adequate description of the experimental conversions, as shown as the solid line in Fig. 6.

Actually, as stressed above, the specificity of ceria promoted Pt/CeO₂/Al₂O₃ is to offer a second route for oxygen activation over the ceria surface in parallel to the main route of oxygen activation on the Pt phase. It can be formalized by adding the following elementary steps to the above mechanism restricted to the Pt phase:

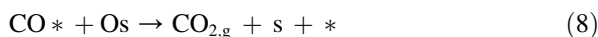

For the present study, due to a lack of kinetic data, the simplified mechanism based on steps (1)–(4) was considered for calculating the CO conversion over the Pt/CeO₂–Al₂O₃ at the reactor outlet. Fig. 6 shows the CO₂ flow as a function of the axial position at 2 different temperatures as calculated by rate Equation (5) (line) and the experimental IR CO₂ bands (points). The adequacy of the chosen simplified model is confirmed for these axial profiles.

One can note here that the above modelling of the CO oxidation along the platelet channels is assuming isothermal conditions, which confirms the excellent heat transfer between the catalytic sites and the metallic substrate.

5. Conclusions

This preliminary study was carried out to demonstrate the potential of an original DRIFTS reactor that was designed for characterizing the “operando”

behaviour of microstructured platelets (element of a stacked microreactor). Its unique ability to map the surface adspecies and gas concentration by means of the DRIFTS/microscope analysis was used for establishing a kinetic model of the CO oxidation able to describe correctly the experimental data corresponding to this mapping.

Comparing the CO conversion changes with temperature between this new tool and a conventional fixed bed DRIFTS reactor brought further evidence of the exceptional heat management of the microstructured reactor (isothermal behaviour without heat transfer limitation) while the fixed bed one was much more adiabatic (leading to quasi autothermal conditions due to major heat transfer limitation).

Thus, despite the drawback of a set-up difficult to handle (combining prototype and commercial units), this tool offers a unique advantage for advanced operando kinetic studies. Further works are in progress to enlarge its field of application.

References

- [1] W. Ehrfeld, V. Hessel, H. Lowe, *Microreactors*, Wiley-VCH, Weinheim, 2000.
- [2] <http://www.atotech.com/>.
- [3] <http://www.velocys.com>.
- [4] A.C. van Veen, Y. Schuurman, C. Mirodatos, *Handbook for Microprocess Engineering*, V. Hessel, A. Renken, J. Schouten, Y. Wang, J.-I. Yoshida (Eds.), *System, Process and Plant Engineering*, vol. C, Wiley VCH, in press.
- [5] M. Bosco, F. Vogel, *Catal. Today* 116 (2006) 348.
- [6] www.Thermoptec.com.
- [7] J. Sirta, S. Phanichphant, F.C. Meunier, *Anal. Chem.* 79 (2007) 3912.
- [8] A. Slagtern, Y. Schuurman, C. Leclercq, X. Verykios, C. Mirodatos, *J. Catal.* 172 (1997) 118.
- [9] F.C. Meunier, A. Goguet, S. Shekhtman, D. Rooney, H. Daly, *Appl. Catal.* 340 (2008) 196.
- [10] A. Stefanescu, A.C. van Veen, C. Mirodatos, J.C. Beziat, E. Duval-Brunel, *Catal. Today* 125 (2007) 16.
- [11] S.P. Teh, Ph.D. Thesis, Lyon University, 24 Nov 2005.
- [12] M.A.J. Campman, Ph.D Thesis, TUE, 1996.
- [13] H. Nibbelke, M.A.J. Campman, J.H.B.J. Hoebink, G.B. Marin, *J. Catal.* 171 (1997) 358.
- [14] G. Germani, P. Alphonse, M. Courty, Y. Schuurman, C. Mirodatos, *Catal. Today* 110 (1–2) (2005) 114.
- [15] A.C. Luntz, M.D. Williams, D.S. Bethune, *J. Chem. Phys.* 89 (1988) 4381.
- [16] R.K. Herz, S.P. Marin, *J. Catal.* 65 (1980) 281.
- [17] V.P. Zhdanov, J. Pavlicek, Z. Knor, *Catal. Rev. Sci. Eng.* 30 (1988) 501.